# Massagekissen

Massage pillow

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# Symbols

E

	Make sure that you read through the operating instructions before first use and observe all instructions.
	Earth-free applied part Body Floating type. The product affords protection against electric shock at permissible leakage current.
	This product bears the recycling symbol according to the EU Guideline 2002/96/EC. Disposal of it via the domestic/residual refuse is not permitted. The contact partner responsible for the disposal is the dealer/manufacturer.
	Manufacturer
	Year of manufacture
	Serial number of the device
	Batch number
	Usable until

# Contents

<b>General informations</b>	<b>3</b>
Device overview	3
Scope of delivery	3
<b>Important notes</b>	<b>4</b>
... on these instructions for use	4
... on the specific function	4
... for liability and warranty	4
Safety instructions	5
<b>Informations</b>	<b>7</b>
Application	7
Massage positions	7
<b>Device description</b>	<b>8</b>
Name and function of the device	8
Technical data and support	8
<b>Care and disposal</b>	<b>9</b>
Cleaning	9
Maintenance	9
Correct Disposal of this product	9

# General informations

E

Dear customer,

many thanks for trusting in us and selecting this top-quality product from

tic Medizintechnik GmbH & Co. KG.

These instructions for use will guide you through the application options of the massage pillow SaneoRELAX.

Please always keep the instruction for use with the device!

---

## Device overview

- 2 programs (massage, massage + heat-function)
- Simple operation
- Velcro fastener for easy mounting

---

## Scope of delivery

1x SaneoRELAX massage pillow

1x Manual

1x Power plug

# Important notes

## ... on these instructions for use

Please read these instructions for use prior to initial use.

These instructions for use form part of the device and must be available to the user at all times. Excerpts are not permitted.

Please carefully read the notes on proper use (S. 4) and the Safety notes (S. 5ff).

## ... on the specific function

The SaneoRELAX is a household appliance for massage with or without heat function.

It may be used only factory authorized accessories.

Any use that is not in the operating instructions or described as not permitted applies as unintended.

Improper use, non-compliance with this document, as well as unauthorized alterations excludes the manufacturer of liability for any resulting property damage and personal injury, and will void the manufacturer's warranty.

## ... for liability and warranty

A warranty period of 24 months applies for the SaneoRELAX massage pillow and the accessories from the day of delivery. During this period the device/accessories can be repaired or replaced free of charge by the manufacturer.

Unintended use, the non-compliance with these operating instructions as well as unauthorised changes carried out of the system preclude the manufacturer's liability for the resulting property damage and personal injury expires the warranty.

The manufacturer reserves the right to make technical changes, which serve the improvement of the therapeutic device, the functionality and the accessories.

# Important notes

## Safety instructions

The manufacturer includes safety checks expressly.

To reduce the risk of electric shock always unplug this product from the electrical outlet immediately after use and before cleaning.

E

To reduce the risk of burns, fire, electric shock, or injury to persons:

- The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or it is dropped into water. Return the product to service center for examination and repair.
- Do not carry this product by its supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never use the product with the air openings blocked. Keep the air openings free of lint, hair and such like.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is administered.

## Important notes

E

- To disconnect, switch off all controls and remove plug from outlet.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not misuse the product.
- Use heated surfaces carefully, they may cause serious burn. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The device may only be used in the home environment.

# Informations

## Application

- Connect the cable's end of the massage pillow with the power plug and put the plug into the fitting plug socket.
- Press the Power Switch once.  
The product will start its kneading massage and heat therapy for 20 minutes.
- Press the Power Switch again.  
This time, the product will stop its heat therapy while the kneading massage continues.
- Press the Power Switch and hold for 3 seconds to turn the product off.
- After use:  
When the product is not in use, it should be unplugged.

## Massage positions

Neck

Abdomen

Calf

Lower Back


Upper Back

Thigh

# Device description

## Name and function of the device
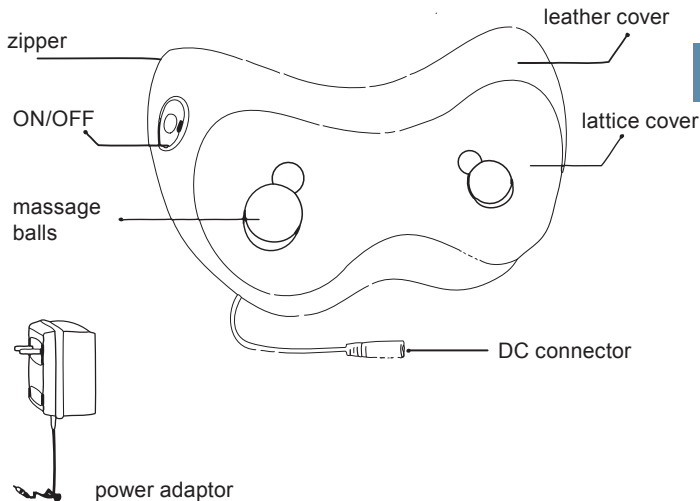
## Technical data and support

Name: SaneoRELAX Massage pillow

Type: ti0066

Adapter Rated Voltage: DC 12V 2.5A

Produkt Rated Power: 30W

tic Medizintechnik GmbH & Co. KG

Endelner Feld 9

46286 Dorsten, Deutschland

Fon: +49 (0) 23 69 . 20 85 - 0

Fax: +49 (0) 23 69 . 20 85 - 10

Web: [www.saneo-line.com](http://www.saneo-line.com)

eMail: [info@saneo-line.com](mailto:info@saneo-line.com)

# Care and disposal

## Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket.

Please use dry cloth to clean the product.

E

### Attention!

Please make sure that do not use benzine or thinner and other aggressive detergents to clean the product as those may damage or let pale the product.

## Maintenance

Make sure that the switch is in the off location and the power plug has been pulled out from the socket.

Keep product clean and far from the wet environment.

## Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

**tic** Medizintechnik GmbH & Co. KG

Endelner Feld 9 · 46286 Dorsten  
Fon: +49 (0) 23 69 · 20 85 - 0  
Fax: +49 (0) 23 69 · 20 85 - 10  
eMail: [info@ticmed.de](mailto:info@ticmed.de)  
Web: [www.ticmed.de](http://www.ticmed.de)

07/2011 6/rt JH