Sport- & Massagerolle

Sports & massage roll

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Contents

General informations	13
Device overview	13
Scope of delivery	13
Important notes	14
... on these instructions for use	14
... on the specific function	14
... for liability and warranty	14
Safety instructions	15
Application	17
Device description	18
Name and function of the device	18
Technical data and support	18
Care and disposal	19
Cleaning	19
Maintenance	19
Battery	19
Correct Disposal of this product	19
Exercise	20

General informations

Dear customer,

many thanks for trusting in us and selecting this top-quality product from

tic Medizintechnik GmbH & Co. KG.

These instructions for use will guide you through the application options of the SaneoROLL+ sports and massage roll.

Please always keep the instruction for use with the device!

E

Device overview

- 4 adjustable intensity levels with LED display
- easy to use with One-Button-Control
- ergonomic design
- hygienically wipeable surface
- time control for a safe application

Scope of delivery

1x SaneoROLL+ sports and massage roll
1x Manual
1x Power plug

Important notes

E

... on these instructions for use

Please read these instructions for use prior to initial use. These instructions for use form part of the device and must be available to the user at all times. Excerpts are not permitted. Please carefully read the notes on proper use and the Safety notes on the following sites.

... on the specific function

The SaneoROLL+ is a household appliance for sports and massage with four intensity levels of depth vibration. The SaneoROLL+ is designed for the application of a pressure and vibration massage by weight to solve muscle tension and muscle regeneration in adults.

It may be used only factory authorized accessories. Any use that is not in the operating instructions or described as not permitted applies as unintended.

... for liability and warranty

A warranty period of 24 months applies for the SaneoROLL+ sports and massage roll and the accessories from the day of delivery. During this period the device/accessories can be repaired or replaced free of charge by the manufacturer.

Unintended use, the non-compliance with these operating instructions as well as unauthorised changes carried out of the system preclude the manufacturer's liability for the resulting property damage and personal injury expires the warranty.

The manufacturer reserves the right to make technical changes, which serve the improvement of the therapeutic device, the functionality and the accessories.

Important notes

E

Safety instructions

- To reduce the danger of an electric shock, disconnect this product from the power supply immediately after charging the battery and before cleaning. During the charging session the power supply has to be adapted in a way it is freely available.
- Before connecting the system to the power supply please check that the specified supply voltage on the name plate matches with the supply voltage in your socket.
- Only charge the system if it is switched off and also only put the power supply into the socket in off mode.
- The product should never be unattended if it is connected with electrical supply. Put the power supply out of the socket if it is not in use.
- Do not use it under a blanket or pillow. An excessive heating up can cause brand and/or electric shock or evoke injuries of persons.
- A surveillance of third parties is necessary if the system is used near to children, invalids or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Never use this product if it has been dropped, has a damaged cable or plug, shows malfunction or had contact to liquids. Send the product to the service center for checking and repairing.
- Do not use the cable as carry handle and keep it away from hot surfaces.
- Do not use the product if the plug openings are blocked. Keep the plug openings free of dust, hairs and the like.
- Do never put objects in the openings.
- Do not use outdoors. Only use in closed rooms or domestic settings.
- Do not operate the system in surroundings where aerosol (spray) products are used or oxygen is provided.
- This device can also be used by persons with limited physical, sensory or mental capacities or lack of experience and knowledge, as well as by children if they are supervised or instructed in the safe use of the device and understood the risks. Children are not allowed to play with the device. Cleaning and maintenance by the user are not to be done by children without supervision.

Important notes

E

- The system is only to be used with the power supply unit which was offered with the device or which is dealer approved.
- The use of the SaneoROLL+ should be omitted or only be done guided by a doctor/therapist if there are problems with the musculoskeletal system or organs.
- In case of pregnancy, implants, pacemakers, circulatory problems, wounds, contusions, inflammations or varicose veins a consultation of the doctor is absolutely necessary before a possible use.
- Do never catch at the device if it has fallen into water. Immediately pull the power plug/power supply unit out of the socket in case it is connected.
- Do not use the SaneoROLL+ as a support or replacement for medical treatments. Symptoms and chronic suffers could get worse in certain circumstances.
- The SaneoROLL+ is not intended for commercial purposes or medical treatments. In case of doubts concerning health please always consult your attending doctor before use.
- Do not use the SaneoROLL+ in damp rooms (e.g. while showering or having a bath)
- Only use the SaneoROLL+ in correct positions like explained on the training card.
- Do not use the device in bed. Do not use the SaneoROLL+ before falling asleep.
- Always switch off the SaneoROLL+ with the function button immediately after use.
- Take care that the cable is never bended, clamped or twisted.
- Avoid the contact of the SaneoROLL+ with sharp and keen objects.
- The vibration massage should be switched off while using the SaneoROLL+ in the head, neck and upper back area.

Important notes


Application

- By a single press of the function button you switch on the SaneoROLL+ and adjust the first vibration level. By a second press you switch to the next higher vibration level. After reaching the fourth intensity level, the device switches off after a further press of the button.
- If the SaneoROLL+ is switched on, it can be switched off by a approx. 2 second pressing of the function button.
- Use the SaneoROLL+ in consideration of the safety instructions according to the training card.
- To charge the battery of the SaneoROLL+, initially switch off the device and then connect it with the power plug. During charging, the status LEDs light up blue. The SaneoROLL+ is fully charged when all four status LEDs are blue.
- The maximum operating period of the SaneoROLL+ for a treatment with vibration is 15 minutes. After that the SaneoROLL+ switches off automatically.

Device description

Name and function of the device

E

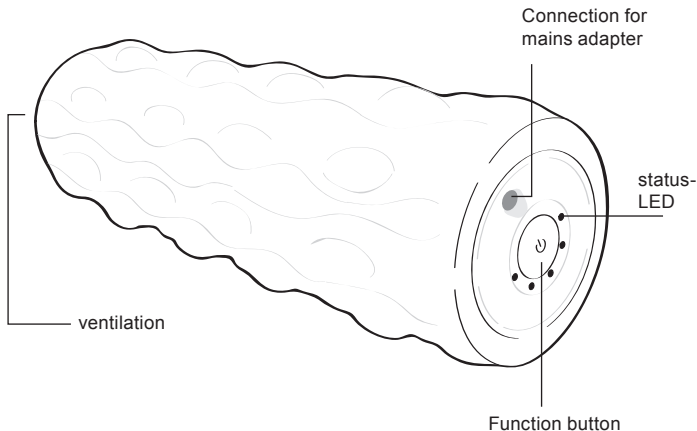Technical data and support

Name: SaneoROLL+ sports and massage roll

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Care and disposal

Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket. Please use dry cloth to clean the product.

Attention!

Please make sure that do not use benzine or thinner and other aggressive detergents to clean the product as those may damage or let pale the product.

Maintenance

Make sure that the switch is in the off location and the power plug has been pulled out from the socket.

Keep product clean and far from the wet environment.

Battery

Don't attempt to replace the SaneoROLL+ battery yourself – you may damage the battery, which could cause overheating, fire, and injury. The lithium-ion battery in your SaneoROLL+ should be serviced or recycled by the manufacturer or an authorized service provider, and must be recycled or disposed of separately from household waste. Dispose of batteries according to your local environmental laws and guidelines.

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Always take care of a correct carrying out of the exercises described here. Use a soft and non-slip mat like e.g. a yoga mat.

Only use the vibration with the indicated exercises and adjust the vibration level to your own feeling.

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Always avoid the use of the SaneoROLL+ in combination with a vibration level in the neck area!

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If this symbol is displayed, perform the exercise without vibration of the SaneoROLL+.

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Place the SaneoROLL+ longitudinally between your shoulder blades. Open your arms at right angles and close them back in front of your chest.

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Position the SaneoROLL+ across your back. Slowly roll up and down on the SaneoROLL+.

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
Lie on the side with your hands behind your head and position the SaneoROLL+ about six inches under your armpit. Keep your lower leg straight and set the other behind it to support your body. Now roll up and down slightly.

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
Put your left hand behind you for support and position the SaneoROLL+ under your left glute. Cross your left leg over your right and then slowly roll back and forth. After about 60 minutes, repeat the exercise on the other side.

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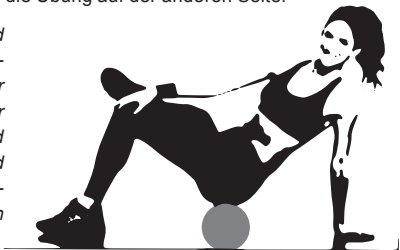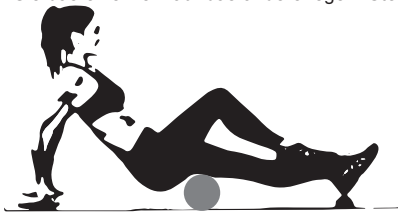Start with the SaneoROLL+ in the middle of your thigh while crossing your opposite leg over the other. Roll down until just before your knee. Now roll up until just your upper muscle insertion of your thigh.

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
Position the SaneoROLL+ on the outside of your right thigh. Place your left leg in front of your right leg and prop yourself on your right hand. Now, slowly roll up and down for approx. one minute. Afterwards change the legs/side and repeat the exercise.

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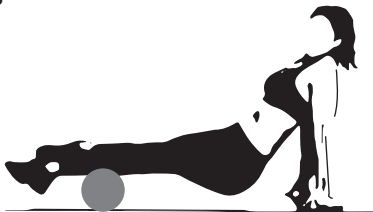Put the SaneoROLL+ in the middle of your left calf muscle. Cross your right leg over your left one. Support yourself now with both hands up. Roll slowly on the SaneoROLL+ back and forth. After 60 seconds repeat the exercise with the other leg.

09/2017/R1 SAN