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Neck and shoulder massager

Saneo+
SENSE

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
Symbols

E

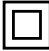	Make sure that you read through the operating instructions before first use and observe all instructions.
	This product bears the recycling symbol according to the EU Guideline 2002/96/EC. Disposal of it via the domestic/residual refuse is not permitted. The contact partner responsible for the disposal is the dealer/manufacturer.
	Manufacturer
	Year of manufacture
	Batch number
	Usable until
	Protection class II - The power supply is subject to protection class II
	Protection class III - SaneoSENSE++ is subject to protection class III

Contents

General informations	3
Device overview	3
Scope of delivery	3
Important notes	4
... on these instructions for use	4
... on the specific function	4
... for liability and warranty	4
Safety instructions	5
Informations	7
Charging the SaneoSENSE+	7
Application	7
Massage positions	7
Device description	8
Name and function of the device	8
Technical data and support	8
Care and disposal	9
Cleaning	9
Maintenance	9
Correct Disposal of this product	9

General informations

E

Dear customer,

many thanks for trusting in us and selecting this top-quality product from

tic Medizintechnik GmbH & Co. KG.

These instructions for use will guide you through the application options of the battery powered SaneoSENSE+ neck and shoulder massager.

Please always keep the instruction for use with the device!

Device overview

- 2 programs (massage, massage + heat-function)
 - 3 Functions (direction change, speed, heat)
 - 3 Speed levels
 - Time control (approx. 15 min.)
 - Mobile use by battery operation
 - Simple operation
-

Scope of delivery

1x SaneoSENSE+ neck and shoulder massager

1x Manual

1x AC adapter

Important notes

... on these instructions for use

Please read these instructions for use prior to initial use.

These instructions for use form part of the device and must be available to the user at all times. Excerpts are not permitted.

Please carefully read the notes on proper use (S. 4) and the Safety notes (S. 5ff).

... on the specific function

The SaneoSENSE+ is a household appliance for massage with or without heat function.

It may be used only factory authorized accessories.

Any use that is not in the operating instructions or described as not permitted applies as unintended.

Improper use, non-compliance with this document, as well as unauthorized alterations excludes the manufacturer of liability for any resulting property damage and personal injury, and will void the manufacturer's warranty.

... for liability and warranty

A warranty period of 24 months applies for the SaneoSENSE+ neck and shoulder massager and the accessories from the day of delivery. During this period the device/accessories can be repaired or replaced free of charge by the manufacturer.

Unintended use, the non-compliance with these operating instructions as well as unauthorised changes carried out of the system preclude the manufacturer's liability for the resulting property damage and personal injury expires the warranty.

The manufacturer reserves the right to make technical changes, which serve the improvement of the therapeutic device, the functionality and the accessories.

Important notes

Safety instructions

The manufacturer includes safety checks expressly.

E To reduce the risk of electric shock, disconnect this product from the power supply immediately after charging the battery and before cleaning. During the charging session the power supply has to be adapted in a way it is freely available.

To reduce the risk of burns, fire, electric shock, or injury to persons:

- The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or it is dropped into water. Return the product to service center for examination and repair.
- Do not carry this product by its supply cord or use cord as a handle.
- Use this product only with the supplied mains adapter.
- Keep the cord away from heated surfaces.
- Never use the product with the air openings blocked. Keep the air openings free of lint, hair and such like.
- Never drop or insert any object into any opening.
- Do not use outdoors.

Important notes

E

- Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- To disconnect the device from the mains adapter, first remove the plug from the socket and then disconnect the device from the mains adapter.
- Operation of the device is only possible if it is disconnected from the mains adapter and the battery is charged.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not misuse the product.
- Use heated surfaces carefully, they may cause serious burn. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The device may only be used in the home environment.
- The device must never be used in or under water.

Informations

Loading the SaneoSENSE +

- Connect the cable end of the SaneoSENSE + to the mains adapter and plug the mains adapter into a suitable socket. The LED lights alternately.
- The SaneoSENSE + battery is fully charged when all the LEDs light up blue. This can take up to three hours.
- After charging, unplug the AC adapter from the AC outlet, and then disconnect the AC adapter from the AC adapter.

Application

- Press the Power Switch once. The product will start its kneading massage and heat therapy for 15 minutes.
- Press the heat button. This time, the product will stop its heat therapy while the kneading massage continues.
- By pressing the button of direction change, the massage balls change their direction of rotation.
- Press the speed button to change the speed of rotation. There are three speed levels which can be set by repeated pressing of the speed button.
- Press the Power Switch and hold for 3 seconds to turn the product off.

Massage positions

Neck

Shoulder


Arms & Legs

Back

Device description

Name and function of the device

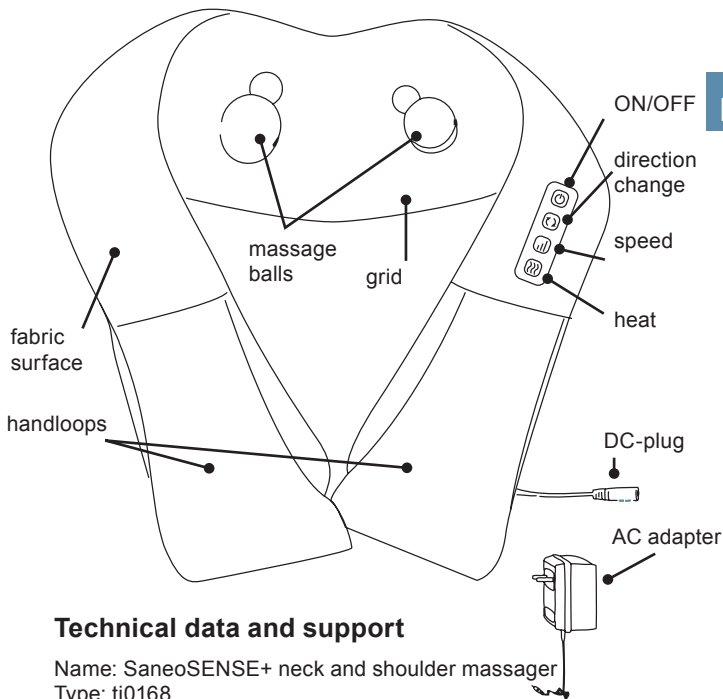Technical data and support

Name: SaneoSENSE+ neck and shoulder massager

Type: ti0168

Adapter Rated Voltage: DC 12.6V/2.2A

Battery Nominal Voltage: 10.8V / 2.2A

Battery Charging Voltage: 12.6V

tic Medizintechnik GmbH & Co. KG

Endelner Feld 9

46286 Dorsten, Deutschland

Fon: +49 (0) 23 69 . 20 85 - 0

Fax: +49 (0) 23 69 . 20 85 - 10

Web: www.saneo-line.com

eMail: info@saneo-line.com

Care and disposal

Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket.

Please use dry cloth to clean the product.

E

Attention!

Please make sure that do not use benzine or thinner and other aggressive detergents to clean the product as those may damage or let pale the product.

Maintenance

Make sure that the switch is in the off location and the power plug has been pulled out from the socket.

Keep product clean and far from the wet environment.

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

tic Medizintechnik GmbH & Co. KG

Endelner Feld 9 · 46286 Dorsten
Fon: +49 (0) 23 69 . 20 85 - 0
Fax: +49 (0) 23 69 . 20 85 - 10
eMail: info@ticmed.de
Web: www.ticmed.de

09/2017/R15AN